

How to activate LIFEtool Solutions Windows applications online?




Prerequisites:


- You have acquired an activation key for the application.
- The application is installed on the computer on which it is to be activated.
- There is an active internet connection on the corresponding computer. If this is not the case, the application must be activated offline.

Procedure:


1. Start the application that is to be activated. As long as the application has not been activated, the licence assistant opens automatically on every start. The assistant informs you how long the application can be used without activation and allows you to activate the application.

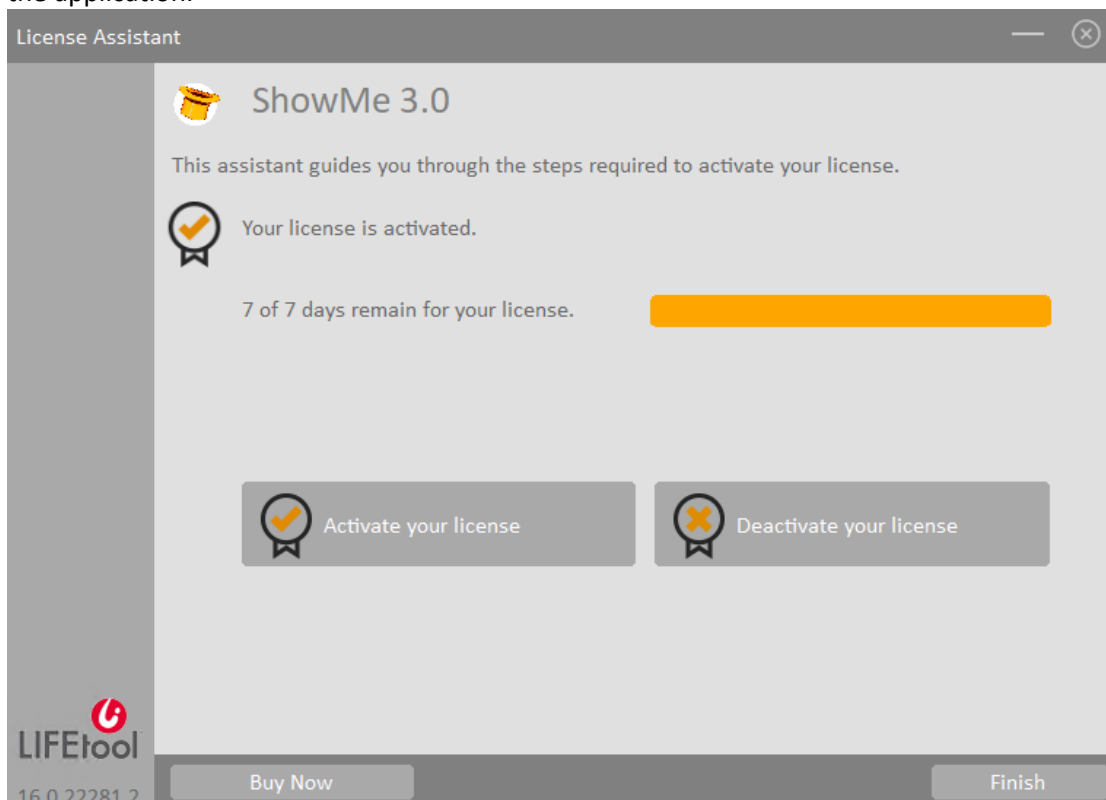


2. To do this, click on "Activate your licence".

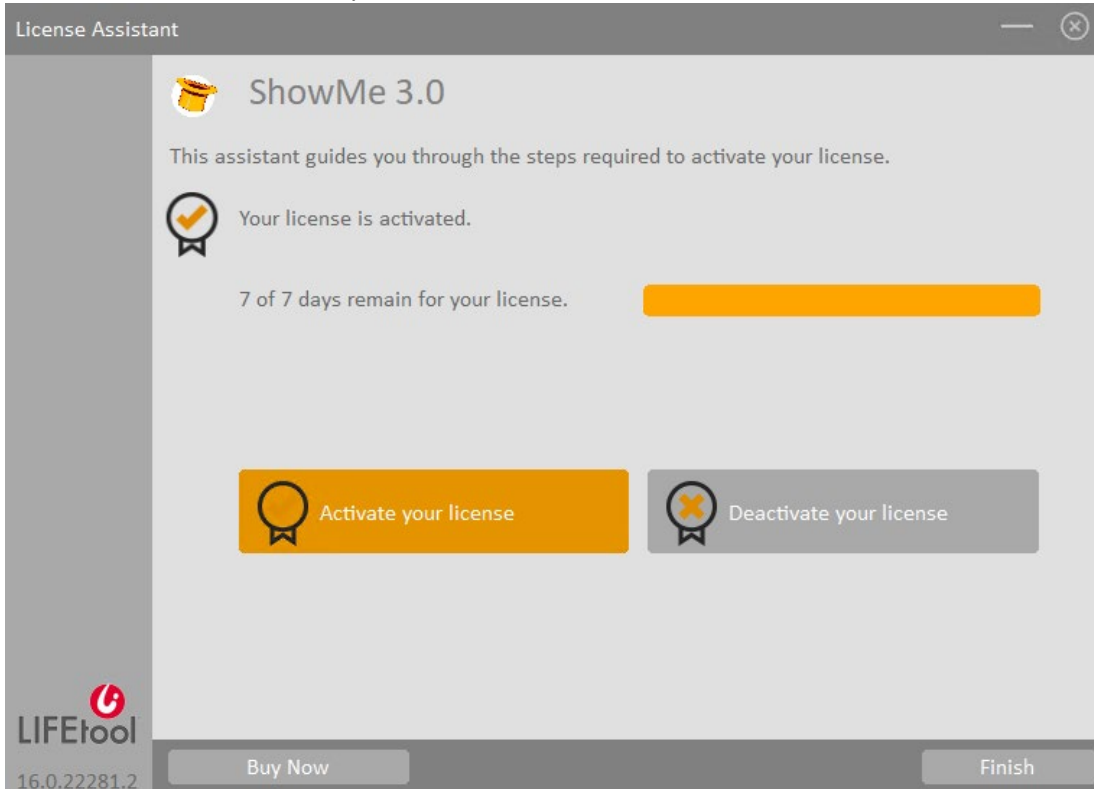


3. In the next window, click "Activate Online".

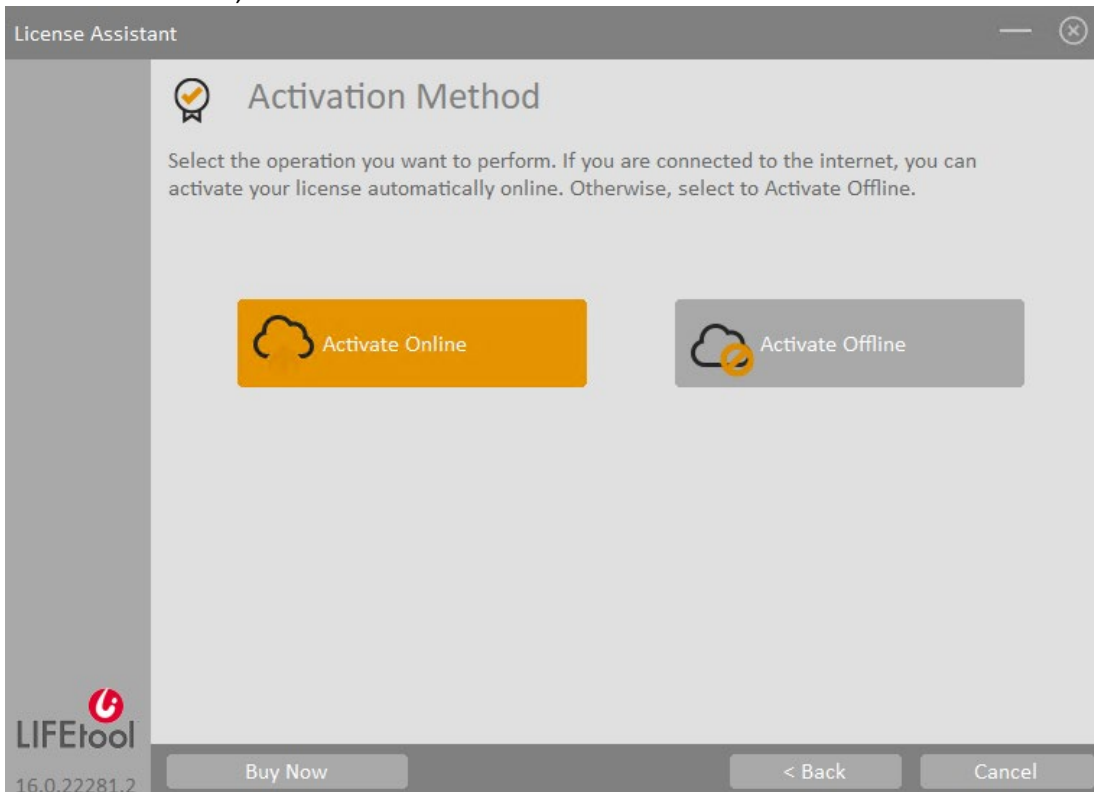

4. Finally, enter the activation key you have received for the application in the "Activation Key" field and click "Activate".

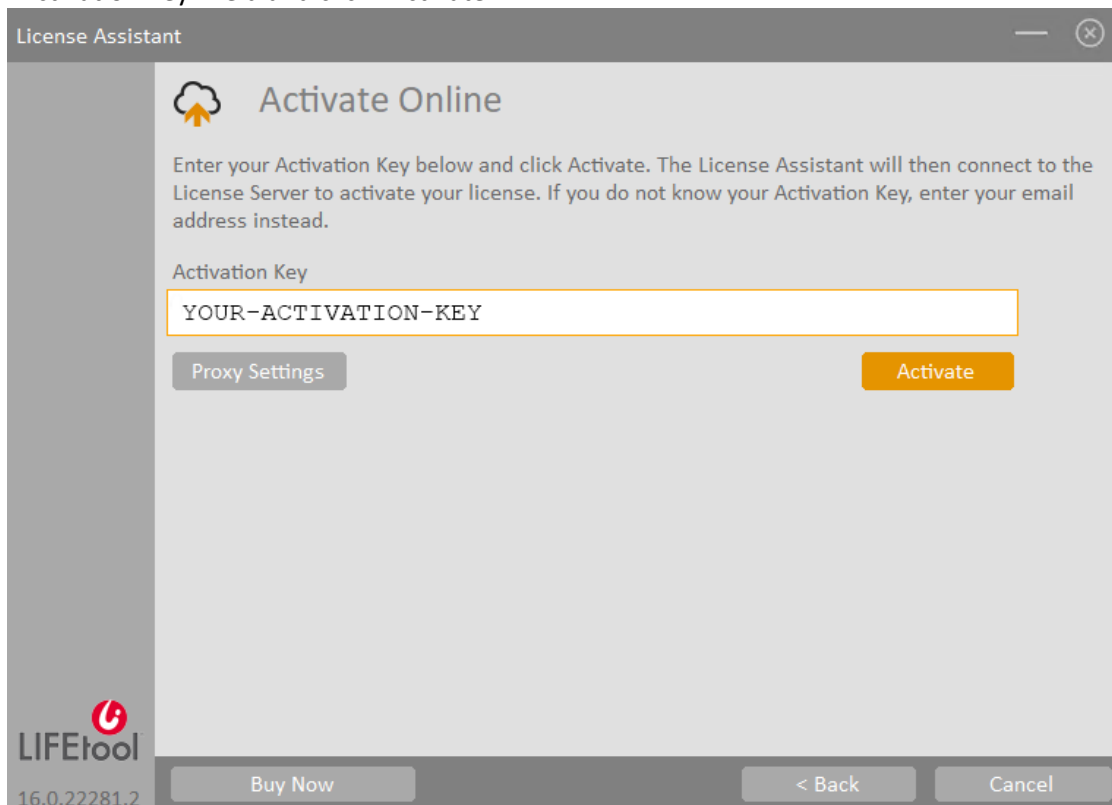

5. Your application is now activated and can be used by clicking "Finish".

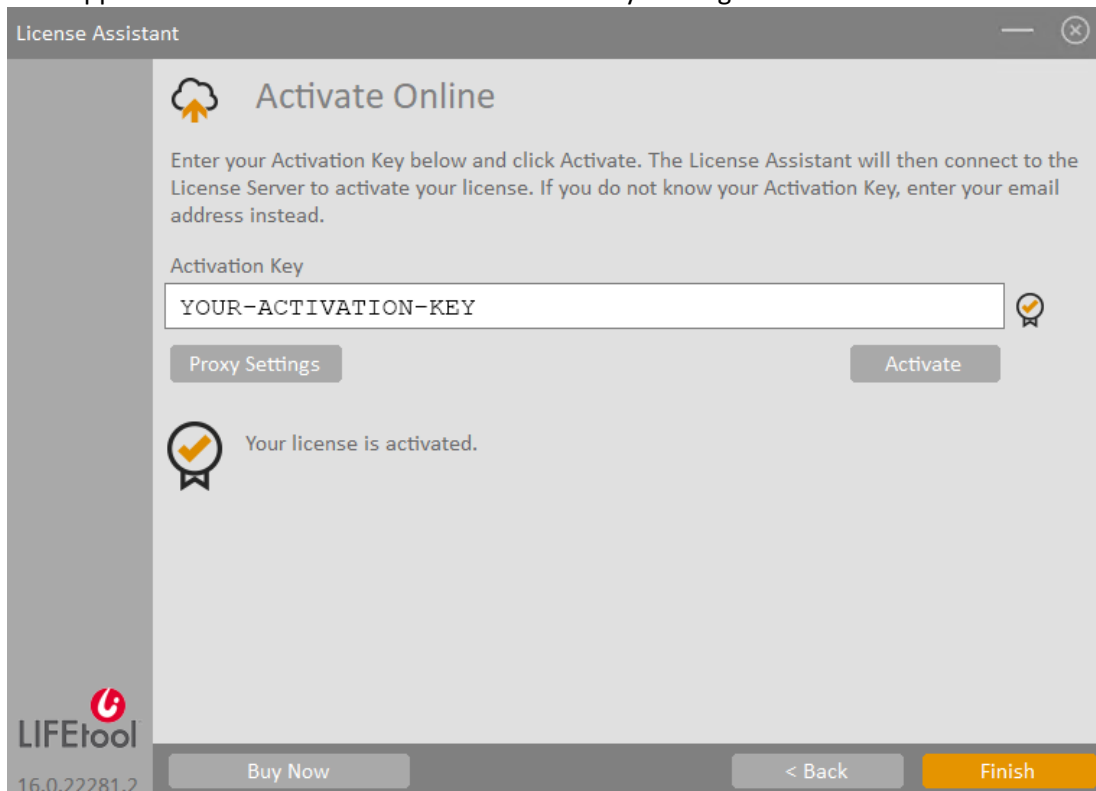

How to activate LIFEtool Solutions Windows applications offline?




If the computer on which the application is to be activated does not have access to the internet or if a firewall blocks the application's access to the internet, it is still possible to activate the application properly.

Prerequisites:


- You have acquired an activation key for the application.
- The application is installed on the computer on which it is to be activated.
- You have a device that is connected to the internet.

Procedure:


1. Start the application that is to be activated. As long as the application has not been activated, the licence assistant opens automatically on every start. The assistant informs you how long the application can be used without activation and allows you to activate the application.

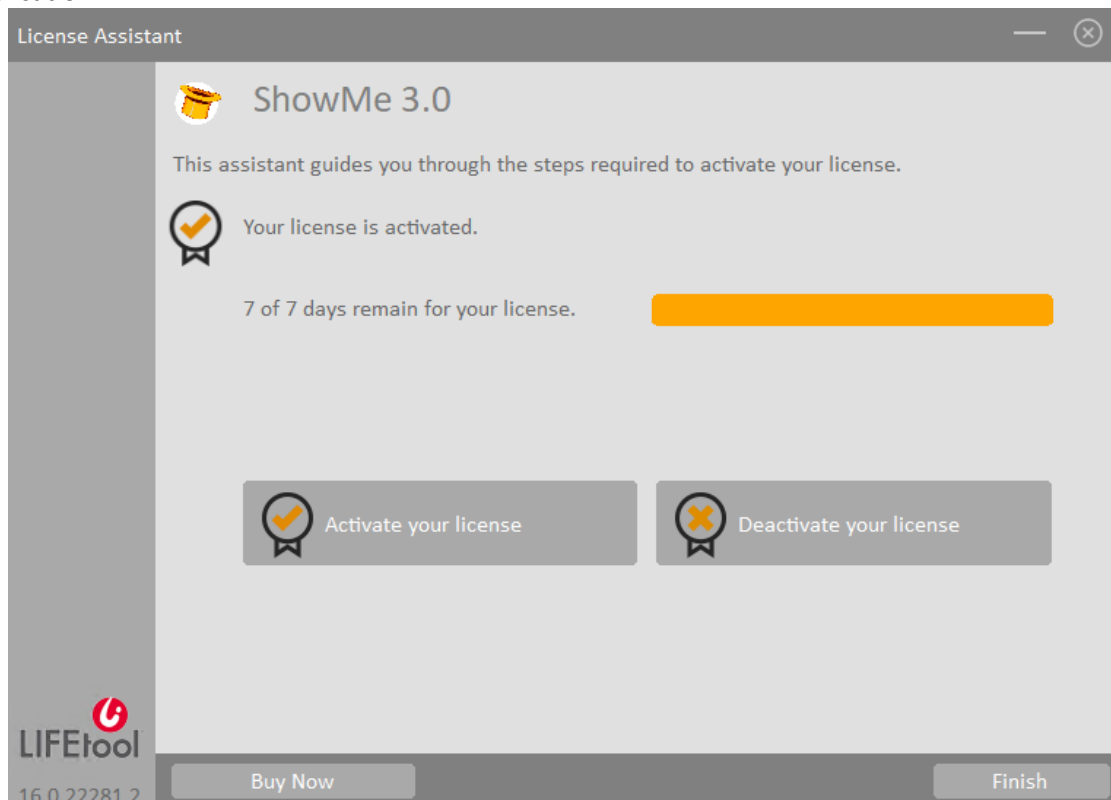


2. To do this, click on "Activate your licence".

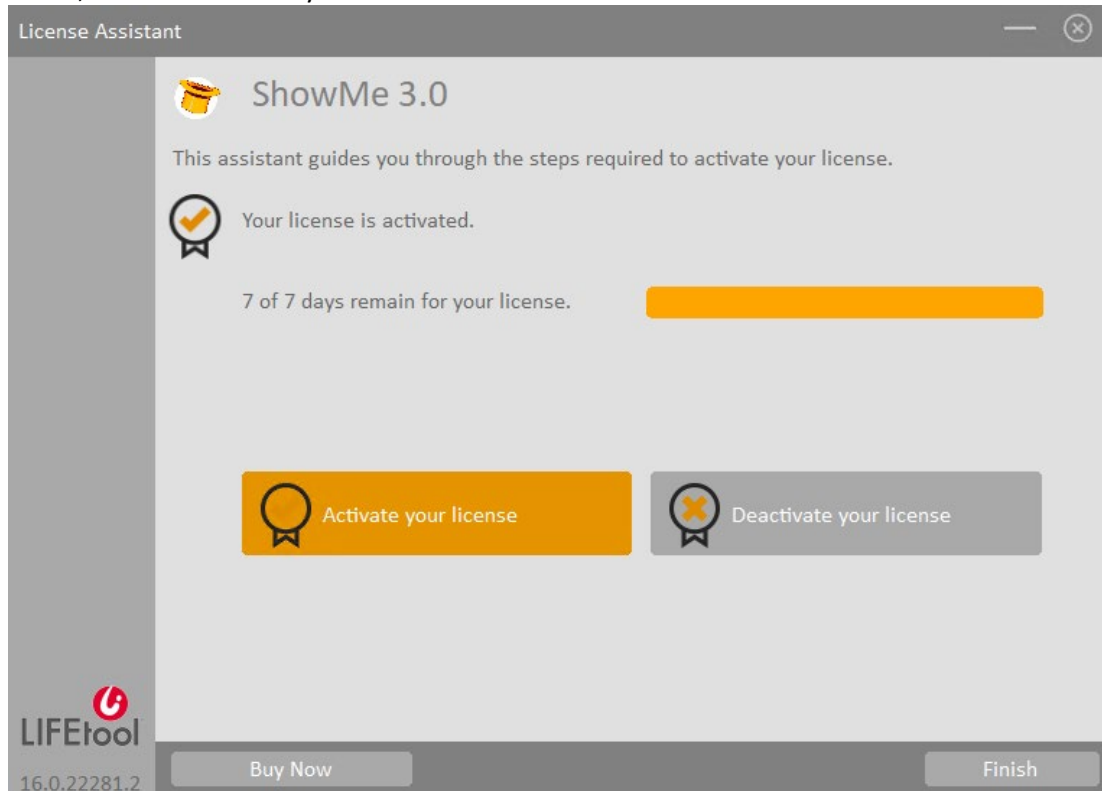


3. In the next window, click "Activate Offline".

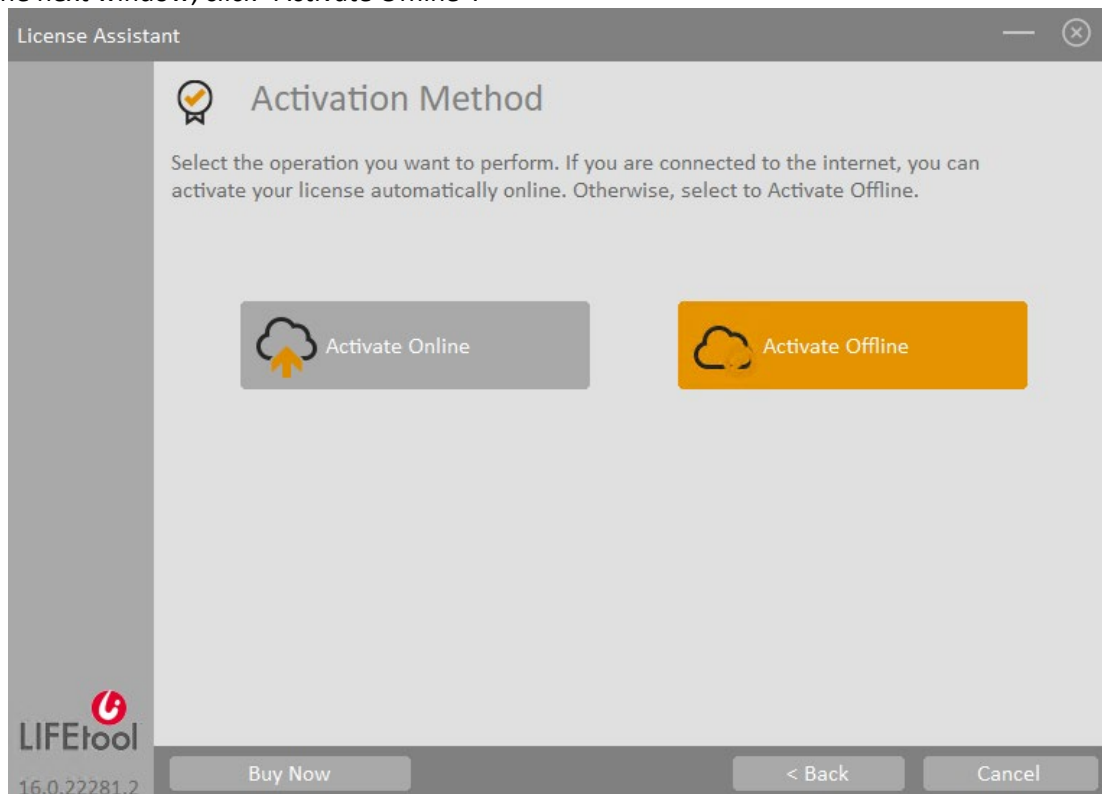

4. Make a note of the text from the "Computer Identifier" field. This will be needed in one of the next steps. If the computer identification is not displayed here, it can be read out from the Windows settings. To do this, open the Windows settings and enter the term "Device name" in the search field and click on the corresponding entry. Make a note of the text that appears under "Device name".

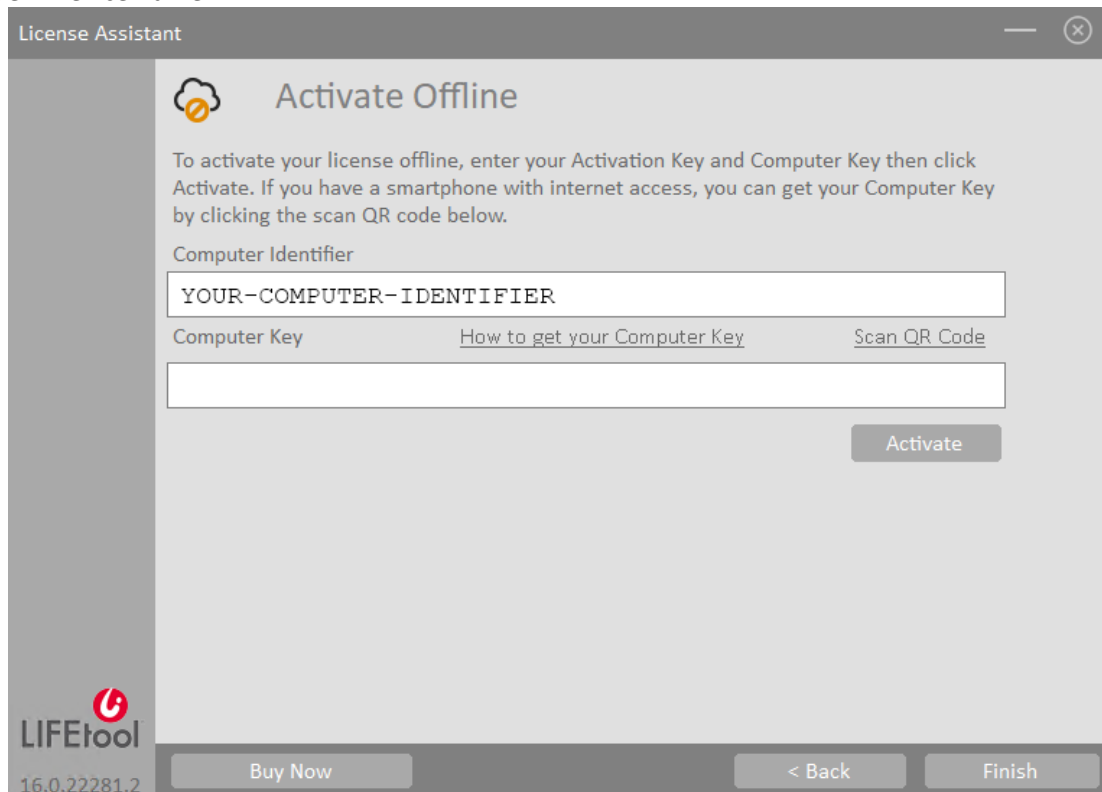

5. Open an internet browser and visit the following page:
<https://qlm3.net/lifetool/qlmcustomersite>

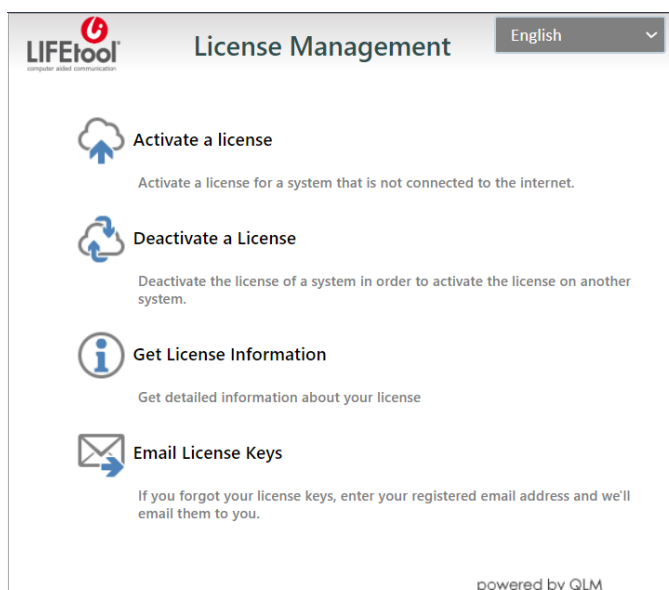

6. Change the language if necessary and click on the field "Activate a licence".

LIFEtool computer aided communication

License Management

English

Activate a license
Activate a license for a system that is not connected to the internet.

Deactivate a License
Deactivate the license of a system in order to activate the license on another system.

Get License Information
Get detailed information about your license

Email License Keys
If you forgot your license keys, enter your registered email address and we'll email them to you.

powered by QLM

7. Now enter your activation key in the first field.

LIFEtool computer aided communication

Offline License Activation

Enter your Activation Key and Computer ID

YOUR-ACTIVATION-KEY


Generate a license file

Computer Key

Activate


powered by QLM

8. For the second field, select the application for which the activation key was created and which is to be activated.

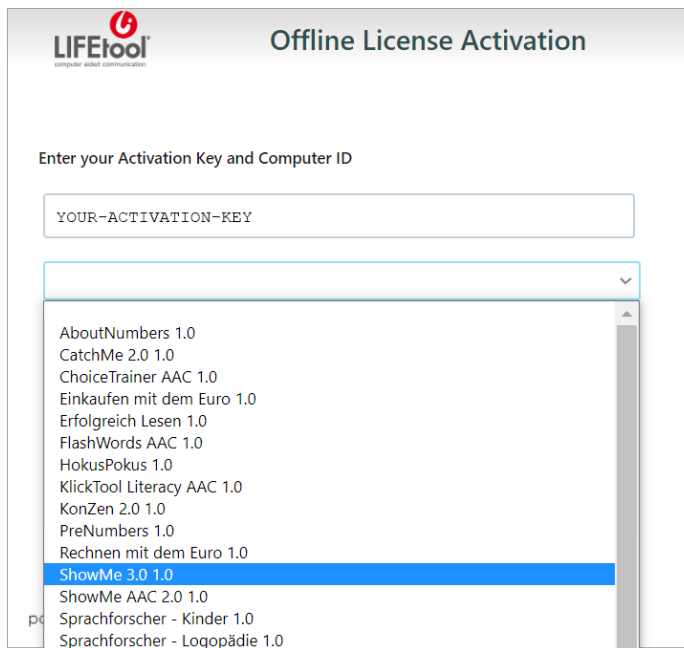

The screenshot shows the 'Offline License Activation' interface. At the top left is the LIFEtool logo with the tagline 'computer aided communication'. The title 'Offline License Activation' is centered at the top. Below the title, the instruction 'Enter your Activation Key and Computer ID' is displayed. There are three input fields: the first contains 'YOUR-ACTIVATION-KEY', the second is a dropdown menu with 'ShowMe 3.0 1.0' selected, and the third contains 'YOUR-COMPUTER-IDENTIFIER'. A list of application names is shown in the dropdown menu, including 'AboutNumbers 1.0', 'CatchMe 2.0 1.0', 'ChoiceTrainer AAC 1.0', 'Einkaufen mit dem Euro 1.0', 'Erfolgreich Lesen 1.0', 'FlashWords AAC 1.0', 'HokusPokus 1.0', 'KlickTool Literacy AAC 1.0', 'KonZen 2.0 1.0', 'PreNumbers 1.0', 'Rechnen mit dem Euro 1.0', 'ShowMe 3.0 1.0', 'ShowMe AAC 2.0 1.0', 'Sprachforscher - Kinder 1.0', and 'Sprachforscher - Logopädie 1.0'. At the bottom left, it says 'powered by QLM'.

9. In the third field, enter the computer ID you noted earlier and then click on the then click on "activate".

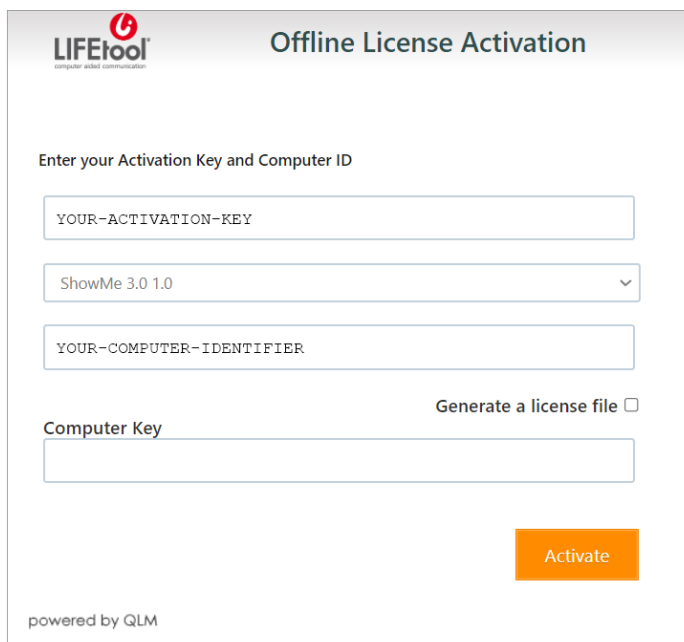

The screenshot shows the 'Offline License Activation' interface. At the top left is the LIFEtool logo with the tagline 'computer aided communication'. The title 'Offline License Activation' is centered at the top. Below the title, the instruction 'Enter your Activation Key and Computer ID' is displayed. There are three input fields: the first contains 'YOUR-ACTIVATION-KEY', the second is a dropdown menu with 'ShowMe 3.0 1.0' selected, and the third contains 'YOUR-COMPUTER-IDENTIFIER'. Below the third field is a checkbox labeled 'Generate a license file' which is unchecked. Below that is a label 'Computer Key' and an empty input field. At the bottom right, there is an orange 'Activate' button. At the bottom left, it says 'powered by QLM'.

10. A computer key will be created, please make a note of it.

The screenshot shows the 'Offline License Activation' window. At the top left is the LIFEtool logo. The title is 'Offline License Activation'. Below the title, it says 'Enter your Activation Key and Computer ID'. There are three input fields: the first contains 'YOUR-ACTIVATION-KEY 961SR-SPKPU', the second is a dropdown menu showing 'ShowMe 3.0 1.0', and the third contains 'YOUR-COMPUTER-IDENTIFIER'. To the right of the third field is a red asterisk. Below these fields is a checkbox labeled 'Generate a license file' which is unchecked. Underneath is a section titled 'Computer Key' with an input field containing 'YOUR-COMPUTER-KEY'. Below this section, there is a message: 'You have successfully activated your license key. Copy the generated Computer Key above and enter it in the corresponding field in the License Wizard.' At the bottom right is an orange 'Activate' button. At the bottom left, it says 'powered by QLM'.

11. Now go back to the licence assistant, enter the noted computer key in the "Computer Key" field and then click "Activate". If the computer identification has not yet been entered, please enter it before clicking "Activate".

The screenshot shows the 'License Assistant' window. The title bar says 'License Assistant'. The main content area has a cloud icon with a checkmark and the title 'Activate Offline'. Below the title, there is a paragraph: 'To activate your license offline, enter your Activation Key and Computer Key then click Activate. If you have a smartphone with internet access, you can get your Computer Key by clicking the scan QR code below.' There are two input fields: the first is labeled 'Computer Identifier' and contains 'YOUR-COMPUTER-IDENTIFIER'; the second is labeled 'Computer Key' and contains 'YOUR-COMPUTER-KEY'. To the right of the 'Computer Key' field are two links: 'How to get your Computer Key' and 'Scan QR Code'. At the bottom right is an orange 'Activate' button. At the bottom left is the LIFEtool logo and the version number '16.0.22281.2'. At the bottom are three buttons: 'Buy Now', '< Back', and 'Cancel'.

12. Your application is now activated and can be started by clicking on "Finish".

